

Shaping the future with you

A collection as diverse as your patients,
Natrelle® is the #1 plastic surgery portfolio.^{1,2,*†}



*Based on surgeon survey data, December 2016 (N = 100).

†Based on US market share data through November 2016 (N = 67).

ART. SCIENCE. *Natrelle*®

Please see inside for Indications and Important Safety Information.

Table of contents

Indications and Important Safety Information	2
<i>Natrelle</i> [®] —the most complete portfolio	8
The <i>Natrelle</i> INSPIRA [®] Collection	13
<i>Natrelle</i> [®] 410 Anatomical Gel Breast Implants	23
<i>Natrelle</i> [®] Saline-Filled Breast Implants	31
<i>Natrelle</i> [®] 133Plus Tissue Expanders	36
<i>Natrelle</i> [®] 133 Tissue Expanders	41
<i>Natrelle</i> [®] tissue expander accessories	46
<i>Natrelle</i> [®] breast implant accessories	47
<i>ConfidencePlus</i> [®] warranty program	49
Additional information	50

**To order breast-shaping products,
please call Customer Care at 1-800-766-0171.**

Lee
Natrelle INSPIRA[®] Style SRM-255
Individual results may vary.

Natrelle® Breast Implants

Important Information

INDICATIONS

Natrelle® Breast Implants are indicated for women for the following:

- **Breast augmentation for women at least 22 years old for silicone-filled implants.**
- **Breast augmentation for women at least 18 years old for saline-filled implants.**

Breast augmentation includes primary breast augmentation to increase breast size, as well as revision surgery to correct or improve the result of a primary breast augmentation surgery.

- **Breast reconstruction.** Breast reconstruction includes primary reconstruction to replace breast tissue that has been removed due to cancer or trauma or that has failed to develop properly due to a severe breast abnormality. Breast reconstruction also includes revision surgery to correct or improve the result of a primary breast reconstruction surgery.

IMPORTANT SAFETY INFORMATION

CONTRAINDICATIONS

Breast implant surgery should not be performed in:

- Women with active infection anywhere in their body.
- Women with existing cancer or precancer of their breast who have not received adequate treatment for those conditions.
- Women who are currently pregnant or nursing.

WARNINGS

- Breast implants are not lifetime devices or necessarily a one-time surgery.
- Avoid damage during surgery: Care should be taken to avoid the use of excessive force and to minimize handling of the implant. Use care when using surgical instruments in proximity with the breast implant. For more information, please see the full Directions for Use.



PRECAUTIONS

Safety and effectiveness have not been established in patients with the following:

- Autoimmune diseases (eg, lupus and scleroderma).
- A compromised immune system (eg, currently receiving immunosuppressive therapy).
- Planned chemotherapy following breast implant placement.
- Planned radiation therapy to the breast following breast implant placement.
- Conditions or medications that interfere with wound healing and blood clotting.
- Reduced blood supply to breast tissue.
- Clinical diagnosis of depression or other mental health disorders, including body dysmorphic disorder and eating disorders. Please discuss any history of mental health disorders prior to surgery. Patients with a diagnosis of depression, or other mental health disorders, should wait until resolution or stabilization of these conditions prior to undergoing breast implantation surgery.

ADVERSE EVENTS

Key adverse events are reoperation, implant removal with or without replacement, implant rupture with silicone-filled implants, implant deflation with saline-filled implants, and capsular contracture Baker Grade III/IV.

Other potential adverse events that may occur with breast implant surgery include: asymmetry, breast pain, breast/skin sensation changes, capsular calcification, delayed wound healing, hematoma, hypertrophic scarring/scarring, implant extrusion, implant malposition, implant palpability/visibility, infection, nipple complications, redness, seroma, swelling, tissue/skin necrosis, wrinkling/rippling.

For more information see the full Directions for Use at www.allergan.com/labeling/usa.htm. To report a problem with Natrelle® Breast Implants, please call Allergan at 1-800-433-8871.

Natrelle® Breast Implants are available by prescription only.

Intraoperative Breast Implant Sizers

Important Information

INDICATIONS

The *Natrelle*® Silicone Sizer and the Allergan Saline Sizer are indicated for single use only for temporary intra-operative insertion in the surgical pocket to evaluate and assist in determining the final breast implant size/volume. The *Natrelle*® Re-sterilizable (410 and Round) Silicone Breast Implant Sizer is used during breast augmentation or reconstruction procedures to assist the surgeon in determining the appropriate size of a breast implant to use.

IMPORTANT SAFETY INFORMATION

CONTRAINDICATIONS

All sizers are contraindicated for use as long-term breast implants or tissue expanders. The *Natrelle*® Silicone Sizer and the Allergan Saline Sizer are contraindicated for multiple patient use or multiple sterilizations.

WARNINGS

Sizers are designed for temporary intra-operative use only and are NOT long-term implants. DO NOT alter, insert or attempt to repair a damaged sizer. DO NOT reuse the *Natrelle*® Silicone Sizer or the Allergan Saline Sizer, which are for single use only. The Silicone Sizers may rupture and release silicone gel. Infection, necrosis, hematoma/seroma and pain may occur following any type of surgery. Minute quantities of silicone gel may diffuse through the elastomer envelope.

PRECAUTIONS

The surgeon must carefully evaluate patient suitability and be knowledgeable about the use of this device. DO NOT expose the sizer to contaminants. Avoid damaging the sizer with surgical instruments (e.g. sharp, blunt or cautery devices). DO NOT attempt to repair damaged products. DO NOT damage the sizer by overhandling, manipulation, or excessive force. Maintain a sterile back-up sizer during surgery.

ADVERSE EVENTS

Adverse events and/or complications may include sepsis, hemorrhage, thrombosis, bleeding, and/or infection.

For more information, please visit www.allergan.com/labeling/usa.htm.

To report a problem, please call Allergan at 1-800-433-8871.

Intraoperative Breast Implant Sizers are available by prescription only.



Natrelle® 133Plus and 133 Tissue Expanders With/Without Suture Tabs and With MAGNA-SITE® Injection Sites

Important Information

INDICATIONS

Natrelle® 133Plus and 133 Tissue Expanders are indicated for:

- Breast reconstruction following mastectomy.
- Treatment of underdeveloped breasts.
- Treatment of soft tissue deformities.

IMPORTANT SAFETY INFORMATION

CONTRAINDICATIONS

Natrelle® 133Plus and 133 Tissue Expanders SHOULD NOT be used in patients:

- Who already have implanted devices that would be affected by a magnetic field (eg, pacemakers, drug infusion devices, artificial sensing devices).
- Whose tissue at the expansion site is determined to be unsuitable.
- Who have an active infection or a residual gross tumor at the expansion site.
- Undergoing adjuvant radiation therapy.
- Whose physiological condition (eg, sensitive over- or underlying anatomy, obesity, smoking, diabetes, autoimmune disease, hypertension, chronic lung or severe cardiovascular disease, or osteogenesis imperfecta) or use of certain drugs (including those that interfere with blood clotting or affect tissue viability) poses an unduly high risk of surgical and/or postoperative complications.
- Who are psychologically unsuitable.

WARNINGS

- DO NOT use Natrelle® 133Plus and 133 Tissue Expanders in patients who already have implanted devices that would be affected by a magnetic field (see Contraindications), because the MAGNA-SITE® integrated injection site contains a strong rare-earth, permanent magnet. Diagnostic testing with Magnetic Resonance Imaging (MRI) is contraindicated in patients with Natrelle® 133Plus and 133 Tissue Expanders in place.
- DO NOT alter the tissue expander or use adulterated fill. Fill only with sterile saline for injection as described in INSTRUCTIONS FOR USE. DO NOT expose to contaminants.
- DO NOT expand if the pressure will compromise wound healing or vasculature of overlying tissue, or beyond patient or tissue tolerance. Stop filling immediately if tissue damage, wound dehiscence, abnormal skin pallor, erythema, edema, pain, or tenderness are observed.
- Active infection anywhere may increase risk of periprosthetic infection. Postoperative infections should be treated aggressively. Unresponsive or necrotizing infection may require premature removal.
- Natrelle® 133Plus and 133 Tissue Expanders are temporary, single-use only devices, and are not to be used for permanent implantation or beyond 6 months. Tissue expansion in breast reconstruction typically requires 4 months to 6 months.
- When using suturing tabs be careful to avoid piercing the shell. Use a new one if damage occurs.

PRECAUTIONS

Active infections may need to be treated and resolved before surgery. Allergan relies on the surgeon to know and follow proper surgical procedures and carefully evaluate patient suitability using standard practice and individual experience. Avoid damage to the tissue expander and use a sterile backup in case of damage. Pay careful attention to tissue tolerance and hemostasis during surgery. Expansion should proceed moderately and never beyond patient or tissue tolerance. Avoid contamination in any postoperative procedure.

ADVERSE REACTIONS

Deflation, tissue damage, infection, extrusion, hematoma/seroma, capsular contracture, premature explantation, displacement, effects on bone, pain, sensation, distortion, inadequate tissue flap, and inflammatory reaction.

**For more information, please visit www.allergan.com/labeling/usa.htm.
To report a problem with Natrelle®, please call Allergan at 1-800-433-8871.**

Natrelle® 133Plus and 133 Tissue Expanders are available by prescription only.



Natrelle® is the #1 plastic surgery portfolio^{1,2,*†}

*Based on surgeon survey data, December 2016 (N = 100).

†Based on US market share data through November 2016 (N = 67).

Artistry becomes reality with the most complete portfolio.³⁻⁵

Countless possibilities.^{6,†}

Inside this catalog, you will find the collection that has it all—the most complete portfolio of cohesive gel breast implants, saline-filled implants, and tissue expanders.³⁻⁵

Our science, research, and development—the substantiating factors in the safety of every one of our products—give us the confidence to know that we are providing you and your patients with a high level of satisfaction.^{7,8}

†Over 700 breast implant options.

Committed to supporting surgeons and patients.

Natrelle® is dedicated to offering you a full range of services for your office and your patients. Programs such as the Natrelle® Gel Rewards Program, Brilliant Distinctions® Consumer Loyalty Program, Natrelle® ConfidencePlus®, and ALLERGAN PARTNER PRIVILEGES® are designed to support your practice and patients every step of the way. And they're only offered from Natrelle® and Allergan.

ART. SCIENCE. *Natrelle*®



The *Natrelle*[®] Cohesive Collection

Natrelle INSPIRA[®] Cohesive

Natrelle[®] 410

Natrelle INSPIRA[®] Soft Touch


Natrelle INSPIRA[®]



Natrelle INSPIRA[®]
ROUND GEL



Natrelle[®] 410
ANATOMICAL GEL

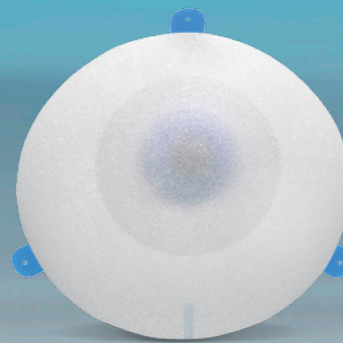

Natrelle[®] 133 PLUS
TISSUE EXPANDER



Natrelle[®] 133
TISSUE EXPANDER

ART. SCIENCE. *Natrelle*[®]

Expansion meets breast shaping

Create a precise pocket for the breast shape
you both envision.

Only *Natrelle*[®] 133Plus Tissue Expanders include the
FOURTE™ Expander Fill System.

- FOURTE™ fills tissue expanders 4X faster than a 21-gauge needle^{9*,†}

Natrelle[®] 133 Tissue Expanders

- 30 years of clinical experience¹⁰

*Clinical significance has not been established.

[†]**Methodology** The FOURTE™ Expander Fill System and 21-gauge needle are attached to 60-cc syringes filled with water. The injector starts the injection at maximum effort while the time is tracked. When the plunger reaches the end of the syringe barrel, the time is stopped and recorded. This is repeated 20 times, recorded, and compared.

Natrelle INSPIRA®
ROUND GEL

Nikki
Natrelle INSPIRA® Style SRF-415
Individual results may vary.





Three cohesivities.
Countless possibilities.^{6,*}

*Over 700 breast implant options.

The Natrelle INSPIRA® Collection

Natrelle® offers the widest range of cohesivities for every breast, patient type, and aesthetic vision—whether it's for reconstruction, revision, or augmentation.¹¹



Three unique silicone gels feature cross-linking technology for 3 distinct cohesivities.




¹**Methodology** Breast implant gel was measured for material properties using the *BTC-2000*™. The *BTC-2000*™ applies a controlled vacuum, or negative pressure, to the gel while measuring the dynamic response of material deformation using a synchronized target laser. The gel for each device (n = 8 per group) was tested at 3 sites, all at or near the apex of the implant (anterior side). From those measurements, the gel cohesivity was calculated and the relative change was determined.

All Natrelle INSPIRA® Round Gel Breast Implants offer:

- Systematic diameter sizing for simplified implant selection
- Smooth or BIOCELL® textured surface
- 5 projection style options
- Matching sizers available

Natrelle INSPIRA® SMOOTH ROUND GEL IMPLANTS

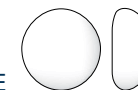

STYLES SCX, SSX, & SRX EXTRA-FULL PROFILE


Smooth Cohesive (SCX)	Smooth Soft Touch (SSX)	Smooth INSPIRA (SRX)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
SCX-200	SSX-200	SRX-200	200	9.50	4.2	MSZ-X205
SCX-225	SSX-225	SRX-225	225	10.00	4.5	MSZ-X225
SCX-255	SSX-255	SRX-255	255	10.25	4.8	MSZ-X255
SCX-285	SSX-285	SRX-285	285	10.50	5.0	MSZ-X285
SCX-310	SSX-310	SRX-310	310	10.75	5.2	MSZ-X310
SCX-340	SSX-340	SRX-340	340	11.00	5.4	MSZ-X340
SCX-375	SSX-375	SRX-375	375	11.25	5.6	MSZ-X375
SCX-400	SSX-400	SRX-400	400	11.50	5.8	MSZ-X400
SCX-420	SSX-420	SRX-420	420	11.75	5.9	MSZ-X420
SCX-445	SSX-445	SRX-445	445	12.00	6.0	MSZ-X445
SCX-470	SSX-470	SRX-470	470	12.25	6.1	MSZ-X470
SCX-495	SSX-495	SRX-495	495	12.50	6.1	MSZ-X495
SCX-525	SSX-525	SRX-525	525	12.75	6.2	MSZ-X525
SCX-545	SSX-545	SRX-545	545	13.00	6.2	MSZ-X545
SCX-560	SSX-560	SRX-560	560	13.25	6.3	MSZ-X560
SCX-580	SSX-580	SRX-580	580	13.50	6.4	MSZ-X580
SCX-615	SSX-615	SRX-615	615	13.75	6.4	MSZ-X615
SCX-650	SSX-650	SRX-650	650	14.00	6.5	MSZ-X650
SCX-700	SSX-700	SRX-700	700	14.25	6.6	MSZ-X700
SCX-750	SSX-750	SRX-750	750	14.50	6.7	MSZ-X750
SCX-800	SSX-800	SRX-800	800	15.50	6.7	MSZ-X800

STYLES SCF, SSF, & SRF FULL PROFILE



Smooth Cohesive (SCF)	Smooth Soft Touch (SSF)	Smooth INSPIRA (SRF)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
SCF-180	SSF-180	SRF-180	180	9.50	4.0	MSZ-F180
SCF-200	SSF-200	SRF-200	200	9.75	4.1	MSZ-F200
SCF-220	SSF-220	SRF-220	220	10.00	4.2	MSZ-F220
SCF-240	SSF-240	SRF-240	240	10.50	4.3	MSZ-F240
SCF-265	SSF-265	SRF-265	265	11.00	4.4	MSZ-F265
SCF-295	SSF-295	SRF-295	295	11.25	4.5	MSZ-F295
SCF-325	SSF-325	SRF-325	325	11.40	4.6	MSZ-F325
SCF-335	SSF-335	SRF-335	335	11.50	4.7	MSZ-F335
SCF-345	SSF-345	SRF-345	345	11.75	4.8	MSZ-F345
SCF-365	SSF-365	SRF-365	365	12.00	4.9	MSZ-F365
SCF-385	SSF-385	SRF-385	385	12.25	5.0	MSZ-F385
SCF-415	SSF-415	SRF-415	415	12.50	5.1	MSZ-F415
SCF-450	SSF-450	SRF-450	450	12.75	5.3	MSZ-F450
SCF-485	SSF-485	SRF-485	485	13.00	5.4	MSZ-F485
SCF-520	SSF-520	SRF-520	520	13.25	5.5	MSZ-F520
SCF-560	SSF-560	SRF-560	560	13.50	5.7	MSZ-F560
SCF-605	SSF-605	SRF-605	605	14.00	5.8	MSZ-F605
SCF-650	SSF-650	SRF-650	650	14.50	5.9	MSZ-F650
SCF-695	SSF-695	SRF-695	695	14.75	6.0	MSZ-F695
SCF-745	SSF-745	SRF-745	745	15.00	6.1	MSZ-F745
SCF-770	SSF-770	SRF-770	770	15.50	6.0	MSZ-F770

Natrelle INSPIRA®
SMOOTH ROUND GEL IMPLANTS

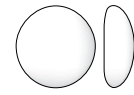


STYLES SCM, SSM, & SRM MODERATE PROFILE

Smooth Cohesive (SCM)	Smooth Soft Touch (SSM)	Smooth INSPIRA (SRM)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
SCM-140	SSM-140	SRM-140	140	9.50	3.0	MSZ-M140
SCM-175	SSM-175	SRM-175	175	10.00	3.3	MSZ-M175
SCM-195	SSM-195	SRM-195	195	10.25	3.4	MSZ-M195
SCM-210	SSM-210	SRM-210	210	10.50	3.5	MSZ-M210
SCM-240	SSM-240	SRM-240	240	11.00	3.7	MSZ-M240
SCM-255	SSM-255	SRM-255	255	11.25	3.8	MSZ-M255
SCM-275	SSM-275	SRM-275	275	11.50	3.9	MSZ-M275
SCM-295	SSM-295	SRM-295	295	11.75	3.9	MSZ-M295
SCM-310	SSM-310	SRM-310	310	12.00	4.0	MSZ-M310
SCM-330	SSM-330	SRM-330	330	12.25	4.1	MSZ-M330
SCM-345	SSM-345	SRM-345	345	12.50	4.2	MSZ-M345
SCM-360	SSM-360	SRM-360	360	12.75	4.2	MSZ-M360
SCM-375	SSM-375	SRM-375	375	13.00	4.3	MSZ-M375
SCM-405	SSM-405	SRM-405	405	13.25	4.4	MSZ-M405
SCM-445	SSM-445	SRM-445	445	13.50	4.5	MSZ-M445
SCM-485	SSM-485	SRM-485	485	14.00	4.6	MSZ-M485
SCM-520	SSM-520	SRM-520	520	14.50	4.8	MSZ-M520
SCM-560	SSM-560	SRM-560	560	14.75	4.9	MSZ-M560
SCM-600	SSM-600	SRM-600	600	15.00	4.9	MSZ-M600
SCM-640	SSM-640	SRM-640	640	15.50	5.0	MSZ-M640
SCM-685	SSM-685	SRM-685	685	16.00	5.2	MSZ-M685
SCM-755	SSM-755	SRM-755	755	16.50	5.0	MSZ-M755

STYLES SCLP, SSLP, & SRLP LOW-PLUS PROFILE

Smooth Cohesive (SCLP)	Smooth Soft Touch (SSLP)	Smooth INSPIRA (SRLP)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
SCLP-125	SSLP-125	SRLP-125	125	9.50	2.5	MSZ-LP125
SCLP-145	SSLP-145	SRLP-145	145	10.00	2.6	MSZ-LP145
SCLP-165	SSLP-165	SRLP-165	165	10.50	2.7	MSZ-LP165
SCLP-190	SSLP-190	SRLP-190	190	11.00	2.8	MSZ-LP190
SCLP-205	SSLP-205	SRLP-205	205	11.25	2.9	MSZ-LP205
SCLP-220	SSLP-220	SRLP-220	220	11.50	3.0	MSZ-LP220
SCLP-235	SSLP-235	SRLP-235	235	11.75	3.1	MSZ-LP235
SCLP-250	SSLP-250	SRLP-250	250	12.00	3.1	MSZ-LP250
SCLP-265	SSLP-265	SRLP-265	265	12.25	3.2	MSZ-LP265
SCLP-280	SSLP-280	SRLP-280	280	12.50	3.3	MSZ-LP280
SCLP-300	SSLP-300	SRLP-300	300	12.75	3.3	MSZ-LP300
SCLP-320	SSLP-320	SRLP-320	320	13.00	3.4	MSZ-LP320
SCLP-340	SSLP-340	SRLP-340	340	13.25	3.5	MSZ-LP340
SCLP-360	SSLP-360	SRLP-360	360	13.50	3.5	MSZ-LP360
SCLP-400	SSLP-400	SRLP-400	400	14.00	3.6	MSZ-LP400
SCLP-440	SSLP-440	SRLP-440	440	14.50	3.8	MSZ-LP440
SCLP-490	SSLP-490	SRLP-490	490	15.00	3.9	MSZ-LP490
SCLP-540	SSLP-540	SRLP-540	540	15.50	4.0	MSZ-LP540
SCLP-590	SSLP-590	SRLP-590	590	16.00	4.2	MSZ-LP590
SCLP-640	SSLP-640	SRLP-640	640	16.50	4.3	MSZ-LP640


Natrelle INSPIRA®
SMOOTH ROUND GEL IMPLANTS

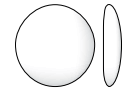

STYLES SCL, SSL, & SRL LOW PROFILE

Smooth Cohesive (SCL)	Smooth Soft Touch (SSL)	Smooth INSPIRA (SRL)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
SCL-110	SSL-110	SRL-110	110	10.00	2.0	MSZ-L110
SCL-125	SSL-125	SRL-125	125	10.25	2.0	MSZ-L125
SCL-140	SSL-140	SRL-140	140	10.50	2.1	MSZ-L140
SCL-170	SSL-170	SRL-170	170	11.00	2.2	MSZ-L170
SCL-200	SSL-200	SRL-200	200	11.50	2.3	MSZ-L200
SCL-230	SSL-230	SRL-230	230	12.00	2.4	MSZ-L230
SCL-260	SSL-260	SRL-260	260	12.50	2.5	MSZ-L260
SCL-290	SSL-290	SRL-290	290	13.00	2.6	MSZ-L290
SCL-320	SSL-320	SRL-320	320	13.50	2.7	MSZ-L320
SCL-350	SSL-350	SRL-350	350	14.00	2.8	MSZ-L350
SCL-380	SSL-380	SRL-380	380	14.50	2.9	MSZ-L380
SCL-410	SSL-410	SRL-410	410	15.00	3.0	MSZ-L410
SCL-460	SSL-460	SRL-460	460	15.50	3.1	MSZ-L460
SCL-510	SSL-510	SRL-510	510	16.00	3.2	MSZ-L510
SCL-560	SSL-560	SRL-560	560	16.50	3.3	MSZ-L560
SCL-610	SSL-610	SRL-610	610	17.00	3.4	MSZ-L610

Natrelle INSPIRA®
TEXTURED ROUND GEL IMPLANTS

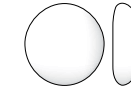


STYLES TCM, TSM, & TRM MODERATE PROFILE

Textured Cohesive (TCM)	Textured Soft Touch (TSM)	Textured INSPIRA (TRM)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
TCM-140	TSM-140	TRM-140	140	9.50	3.0	MSZ-M140
TCM-175	TSM-175	TRM-175	175	10.00	3.3	MSZ-M175
TCM-195	TSM-195	TRM-195	195	10.25	3.4	MSZ-M195
TCM-210	TSM-210	TRM-210	210	10.50	3.5	MSZ-M210
TCM-240	TSM-240	TRM-240	240	11.00	3.7	MSZ-M240
TCM-255	TSM-255	TRM-255	255	11.25	3.8	MSZ-M255
TCM-275	TSM-275	TRM-275	275	11.50	3.9	MSZ-M275
TCM-295	TSM-295	TRM-295	295	11.75	3.9	MSZ-M295
TCM-310	TSM-310	TRM-310	310	12.00	4.0	MSZ-M310
TCM-330	TSM-330	TRM-330	330	12.25	4.1	MSZ-M330
TCM-345	TSM-345	TRM-345	345	12.50	4.2	MSZ-M345
TCM-360	TSM-360	TRM-360	360	12.75	4.2	MSZ-M360
TCM-375	TSM-375	TRM-375	375	13.00	4.3	MSZ-M375
TCM-405	TSM-405	TRM-405	405	13.25	4.4	MSZ-M405
TCM-445	TSM-445	TRM-445	445	13.50	4.5	MSZ-M445
TCM-485	TSM-485	TRM-485	485	14.00	4.6	MSZ-M485
TCM-520	TSM-520	TRM-520	520	14.50	4.8	MSZ-M520
TCM-560	TSM-560	TRM-560	560	14.75	4.9	MSZ-M560
TCM-600	TSM-600	TRM-600	600	15.00	4.9	MSZ-M600
TCM-640	TSM-640	TRM-640	640	15.50	5.0	MSZ-M640
TCM-685	TSM-685	TRM-685	685	16.00	5.2	MSZ-M685
TCM-755	TSM-755	TRM-755	755	16.50	5.2	MSZ-M755

STYLES TCLP, TSLP, & TRLP LOW-PLUS PROFILE

Textured Cohesive (TCLP)	Textured Soft Touch (TSLP)	Textured INSPIRA (TRLP)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
TCLP-125	TSLP-125	TRLP-125	125	9.50	2.5	MSZ-LP125
TCLP-145	TSLP-145	TRLP-145	145	10.00	2.6	MSZ-LP145
TCLP-165	TSLP-165	TRLP-165	165	10.50	2.7	MSZ-LP165
TCLP-190	TSLP-190	TRLP-190	190	11.00	2.8	MSZ-LP190
TCLP-205	TSLP-205	TRLP-205	205	11.25	2.9	MSZ-LP205
TCLP-220	TSLP-220	TRLP-220	220	11.50	3.0	MSZ-LP220
TCLP-235	TSLP-235	TRLP-235	235	11.75	3.1	MSZ-LP235
TCLP-250	TSLP-250	TRLP-250	250	12.00	3.1	MSZ-LP250
TCLP-265	TSLP-265	TRLP-265	265	12.25	3.2	MSZ-LP265
TCLP-280	TSLP-280	TRLP-280	280	12.50	3.3	MSZ-LP280
TCLP-300	TSLP-300	TRLP-300	300	12.75	3.3	MSZ-LP300
TCLP-320	TSLP-320	TRLP-320	320	13.00	3.4	MSZ-LP320
TCLP-340	TSLP-340	TRLP-340	340	13.25	3.5	MSZ-LP340
TCLP-360	TSLP-360	TRLP-360	360	13.50	3.5	MSZ-LP360
TCLP-400	TSLP-400	TRLP-400	400	14.00	3.6	MSZ-LP400
TCLP-440	TSLP-440	TRLP-440	440	14.50	3.8	MSZ-LP440
TCLP-490	TSLP-490	TRLP-490	490	15.00	3.9	MSZ-LP490
TCLP-540	TSLP-540	TRLP-540	540	15.50	4.0	MSZ-LP540
TCLP-590	TSLP-590	TRLP-590	590	16.00	4.2	MSZ-LP590
TCLP-640	TSLP-640	TRLP-640	640	16.50	4.3	MSZ-LP640


Natrelle INSPIRA®
TEXTURED ROUND GEL IMPLANTS

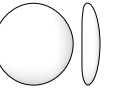

STYLES TCL, TSL, & TRL LOW PROFILE

Textured Cohesive (TCL)	Textured Soft Touch (TSL)	Textured INSPIRA (TRL)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
TCL-110	TSL-110	TRL-110	110	10.00	1.9	MSZ-L110
TCL-125	TSL-125	TRL-125	125	10.50	2.0	MSZ-L125
TCL-140	TSL-140	TRL-140	140	11.00	2.0	MSZ-L140
TCL-170	TSL-170	TRL-170	170	11.50	2.2	MSZ-L170
TCL-200	TSL-200	TRL-200	200	12.00	2.4	MSZ-L200
TCL-230	TSL-230	TRL-230	230	12.50	2.5	MSZ-L230
TCL-260	TSL-260	TRL-260	260	13.00	2.6	MSZ-L260
TCL-290	TSL-290	TRL-290	290	13.50	2.7	MSZ-L290
TCL-320	TSL-320	TRL-320	320	14.00	2.8	MSZ-L320
TCL-350	TSL-350	TRL-350	350	14.25	2.9	MSZ-L350
TCL-380	TSL-380	TRL-380	380	14.50	3.0	MSZ-L380
TCL-410	TSL-410	TRL-410	410	15.00	3.1	MSZ-L410
TCL-460	TSL-460	TRL-460	460	15.50	3.2	MSZ-L460
TCL-510	TSL-510	TRL-510	510	16.00	3.4	MSZ-L510
TCL-560	TSL-560	TRL-560	560	16.50	3.5	MSZ-L560
TCL-610	TSL-610	TRL-610	610	17.00	3.6	MSZ-L610


Natrelle INSPIRA®
TEXTURED ROUND GEL IMPLANTS

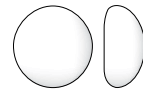

STYLES TCX, TSX, & TRX

EXTRA-FULL PROFILE

Textured Cohesive (TCX)	Textured Soft Touch (TSX)	Textured INSPIRA (TRX)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
TCX-205	TSX-205	TRX-205	205	9.50	4.4	MSZ-X205
TCX-225	TSX-225	TRX-225	225	10.00	4.6	MSZ-X225
TCX-255	TSX-255	TRX-255	255	10.00	4.8	MSZ-X255
TCX-285	TSX-285	TRX-285	285	10.50	5.0	MSZ-X285
TCX-310	TSX-310	TRX-310	310	11.00	5.4	MSZ-X310
TCX-340	TSX-340	TRX-340	340	11.25	5.5	MSZ-X340
TCX-375	TSX-375	TRX-375	375	11.50	5.7	MSZ-X375
TCX-400	TSX-400	TRX-400	400	11.75	5.8	MSZ-X400
TCX-420	TSX-420	TRX-420	420	12.00	5.9	MSZ-X420
TCX-445	TSX-445	TRX-445	445	12.25	6.0	MSZ-X445
TCX-470	TSX-470	TRX-470	470	12.50	6.2	MSZ-X470
TCX-495	TSX-495	TRX-495	495	12.75	6.2	MSZ-X495
TCX-525	TSX-525	TRX-525	525	13.00	6.3	MSZ-X525
TCX-545	TSX-545	TRX-545	545	13.25	6.3	MSZ-X545
TCX-560	TSX-560	TRX-560	560	13.50	6.4	MSZ-X560
TCX-580	TSX-580	TRX-580	580	13.75	6.5	MSZ-X580
TCX-615	TSX-615	TRX-615	615	14.00	6.5	MSZ-X615
TCX-650	TSX-650	TRX-650	650	14.25	6.6	MSZ-X650
TCX-700	TSX-700	TRX-700	700	14.50	6.6	MSZ-X700
TCX-750	TSX-750	TRX-750	750	15.00	6.7	MSZ-X750
TCX-800	TSX-800	TRX-800	800	15.50	6.8	MSZ-X800

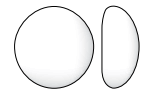

STYLES TCF, TSF, & TRF

FULL PROFILE

Textured Cohesive (TCF)	Textured Soft Touch (TSF)	Textured INSPIRA (TRF)	Volume (cc)	Diameter (cm)	Projection (cm)	Re-sterilizable Sizer
TCF-180	TSF-180	TRF-180	180	9.50	4.0	MSZ-F180
TCF-200	TSF-200	TRF-200	200	9.75	4.1	MSZ-F200
TCF-220	TSF-220	TRF-220	220	10.00	4.2	MSZ-F220
TCF-240	TSF-240	TRF-240	240	10.50	4.3	MSZ-F240
TCF-265	TSF-265	TRF-265	265	11.00	4.5	MSZ-F265
TCF-295	TSF-295	TRF-295	295	11.25	4.6	MSZ-F295
TCF-325	TSF-325	TRF-325	325	11.50	4.7	MSZ-F325
TCF-335	TSF-335	TRF-335	335	11.75	4.8	MSZ-F335
TCF-345	TSF-345	TRF-345	345	12.00	4.9	MSZ-F345
TCF-365	TSF-365	TRF-365	365	12.25	5.0	MSZ-F365
TCF-385	TSF-385	TRF-385	385	12.50	5.1	MSZ-F385
TCF-415	TSF-415	TRF-415	415	12.75	5.2	MSZ-F415
TCF-450	TSF-450	TRF-450	450	13.00	5.3	MSZ-F450
TCF-485	TSF-485	TRF-485	485	13.25	5.4	MSZ-F485
TCF-520	TSF-520	TRF-520	520	13.50	5.5	MSZ-F520
TCF-560	TSF-560	TRF-560	560	14.00	5.7	MSZ-F560
TCF-605	TSF-605	TRF-605	605	14.50	5.8	MSZ-F605
TCF-650	TSF-650	TRF-650	650	14.75	5.9	MSZ-F650
TCF-695	TSF-695	TRF-695	695	15.00	6.0	MSZ-F695
TCF-745	TSF-745	TRF-745	745	15.50	6.1	MSZ-F745
TCF-770	TSF-770	TRF-770	770	16.00	6.2	MSZ-F770

Natrelle® **410**
 ANATOMICAL GEL

Kortnie
 Natrelle® 410 Style FM
 Individual results may vary.





Natrelle® 410 Anatomical Gel Breast Implants

Designed to mirror the curve of a woman's breast.

Natrelle® 410 Anatomical Gel Breast Implants are filled with highly cohesive silicone gel and are designed to maintain their shape over time.


Natrelle® 410 Anatomical Gel Breast Implants have the most high-projection (≥ 6.0 cm) sizes available.³⁻⁵

Natrelle® 410 Anatomical Gel Breast Implants offer:


- Shaped fullness
- Highly cohesive gel
- BIOCELL® textured surface
- 4 projection style options
- Matching sizes available

116 size options³

Offering the highest-projection implant—7.1 cm—in all 3 heights.³⁻⁵

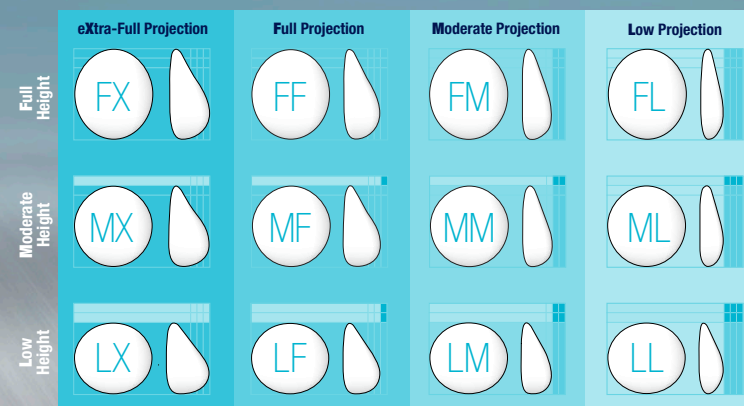

Natrelle® 410 ANATOMICAL GEL IMPLANTS

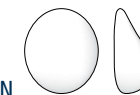

STYLE FX FULL HEIGHT/EXTRA-FULL PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
FX-410185	185	9.5	10.0	4.6	MSZFX185
FX-410215	215	10.0	10.5	4.9	MSZFX215
FX-410245	245	10.5	11.0	5.1	MSZFX245
FX-410280	280	11.0	11.5	5.3	MSZFX280
FX-410315	315	11.5	12.0	5.5	MSZFX315
FX-410360	360	12.0	12.5	5.7	MSZFX360
FX-410410	410	12.5	13.0	6.0	MSZFX410
FX-410450	450	13.0	13.5	6.1	MSZFX450
FX-410495	495	13.5	14.0	6.2	MSZFX495
FX-410560	560	14.0	14.5	6.5	MSZFX560
FX-410615	615	14.5	15.0	6.7	MSZFX615
FX-410690	690	15.0	15.5	7.0	MSZFX690
FX-410775	775	15.5	16.0	7.1	MSZFX775


STYLE MX MODERATE HEIGHT/EXTRA-FULL PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
MX-410165	165	9.5	8.6	4.6	MSZMX165
MX-410195	195	10.0	9.1	4.9	MSZMX195
MX-410225	225	10.5	9.6	5.1	MSZMX225
MX-410255	255	11.0	10.1	5.3	MSZMX255
MX-410290	290	11.5	10.6	5.5	MSZMX290
MX-410325	325	12.0	11.1	5.7	MSZMX325
MX-410370	370	12.5	11.6	6.0	MSZMX370
MX-410410	410	13.0	12.1	6.1	MSZMX410
MX-410445	445	13.5	12.5	6.2	MSZMX445
MX-410520	520	14.0	12.9	6.5	MSZMX520
MX-410550	550	14.5	13.2	6.7	MSZMX550
MX-410620	620	15.0	13.6	7.0	MSZMX620
MX-410685	685	15.5	13.9	7.1	MSZMX685

STYLE LX LOW HEIGHT/EXTRA-FULL PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
LX-410145	145	9.5	7.6	4.6	MSZLX145
LX-410175	175	10.0	8.1	4.9	MSZLX175
LX-410195	195	10.5	8.6	5.1	MSZLX195
LX-410225	225	11.0	9.1	5.3	MSZLX225
LX-410255	255	11.5	9.6	5.5	MSZLX255
LX-410290	290	12.0	10.1	5.7	MSZLX290
LX-410330	330	12.5	10.5	6.0	MSZLX330
LX-410365	365	13.0	10.9	6.1	MSZLX365
LX-410405	405	13.5	11.4	6.2	MSZLX405
LX-410455	455	14.0	11.8	6.5	MSZLX455
LX-410515	515	14.5	12.2	6.7	MSZLX515
LX-410570	570	15.0	12.6	7.0	MSZLX570
LX-410625	625	15.5	13.0	7.1	MSZLX625

Natrelle® 410
ANATOMICAL GEL IMPLANTS

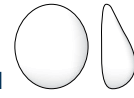

STYLE FF FULL HEIGHT/FULL PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
FF-410185	185	10.0	10.5	4.0	MSZFF185
FF-410220	220	10.5	11.0	4.2	MSZFF220
FF-410255	255	11.0	11.5	4.4	MSZFF255
FF-410290	290	11.5	12.0	4.6	MSZFF290
FF-410335	335	12.0	12.5	4.8	MSZFF335
FF-410375	375	12.5	13.0	5.1	MSZFF375
FF-410425	425	13.0	13.5	5.2	MSZFF425
FF-410475	475	13.5	14.0	5.3	MSZFF475
FF-410535	535	14.0	14.5	5.6	MSZFF535
FF-410595	595	14.5	15.0	5.8	MSZFF595
FF-410655	655	15.0	15.5	6.1	MSZFF655
FF-410740	740	15.5	16.0	6.2	MSZFF740


STYLE MF MODERATE HEIGHT/FULL PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
MF-410140	140	9.5	8.6	3.7	MSZMF140
MF-410165	165	10.0	9.1	4.0	MSZMF165
MF-410195	195	10.5	9.6	4.2	MSZMF195
MF-410225	225	11.0	10.1	4.4	MSZMF225
MF-410255	255	11.5	10.6	4.6	MSZMF255
MF-410295	295	12.0	11.1	4.8	MSZMF295
MF-410335	335	12.5	11.6	5.1	MSZMF335
MF-410375	375	13.0	12.1	5.2	MSZMF375
MF-410420	420	13.5	12.5	5.3	MSZMF420
MF-410470	470	14.0	12.9	5.6	MSZMF470
MF-410525	525	14.5	13.2	5.8	MSZMF525
MF-410580	580	15.0	13.6	6.1	MSZMF580
MF-410640	640	15.5	13.9	6.2	MSZMF640

STYLE LF LOW HEIGHT/FULL PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
LF-410125	125	9.5	7.6	3.7	MSZLF125
LF-410150	150	10.0	8.1	4.0	MSZLF150
LF-410175	175	10.5	8.6	4.2	MSZLF175
LF-410205	205	11.0	9.1	4.4	MSZLF205
LF-410240	240	11.5	9.6	4.6	MSZLF240
LF-410270	270	12.0	10.1	4.8	MSZLF270
LF-410310	310	12.5	10.5	5.1	MSZLF310
LF-410390	390	13.5	11.4	5.3	MSZLF390
LF-410440	440	14.0	11.8	5.6	MSZLF440
LF-410490	490	14.5	12.2	5.8	MSZLF490
LF-410540	540	15.0	12.6	6.1	MSZLF540
LF-410595	595	15.5	13.0	6.2	MSZLF595

Natrelle® 410
ANATOMICAL GEL IMPLANTS

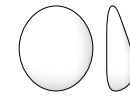

STYLE FM FULL HEIGHT/MODERATE PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
FM-410205	205	10.5	11.0	3.8	MSZFM205
FM-410235	235	11.0	11.5	4.0	MSZFM235
FM-410270	270	11.5	12.0	4.2	MSZFM270
FM-410310	310	12.0	12.5	4.4	MSZFM310
FM-410350	350	12.5	13.0	4.6	MSZFM350
FM-410395	395	13.0	13.5	4.8	MSZFM395
FM-410440	440	13.5	14.0	5.0	MSZFM440
FM-410500	500	14.0	14.5	5.2	MSZFM500
FM-410550	550	14.5	15.0	5.4	MSZFM550
FM-410605	605	15.0	15.5	5.5	MSZFM605
FM-410670	670	15.5	16.0	5.6	MSZFM670


STYLE MM MODERATE HEIGHT/MODERATE PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
MM-410160	160	10.0	9.1	3.6	MSZMM160
MM-410185	185	10.5	9.6	3.8	MSZMM185
MM-410215	215	11.0	10.1	4.0	MSZMM215
MM-410245	245	11.5	10.6	4.2	MSZMM245
MM-410280	280	12.0	11.1	4.4	MSZMM280
MM-410320	320	12.5	11.6	4.6	MSZMM320
MM-410360	360	13.0	12.1	4.8	MSZMM360
MM-410400	400	13.5	12.5	5.0	MSZMM400
MM-410450	450	14.0	12.9	5.2	MSZMM450

STYLE LM LOW HEIGHT/MODERATE PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
LM-410140	140	10.0	8.1	3.6	MSZLM140
LM-410190	190	11.0	9.1	4.0	MSZLM190
LM-410220	220	11.5	9.6	4.2	MSZLM220
LM-410250	250	12.0	10.1	4.4	MSZLM250
LM-410320	320	13.0	10.9	4.8	MSZLM320

Natrelle® 410
ANATOMICAL GEL IMPLANTS

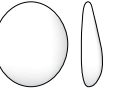

STYLE FL FULL HEIGHT/LOW PROJECTION


Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
FL-410140	140	10.0	10.5	2.9	MSZFL140
FL-410190	190	11.0	11.5	3.1	MSZFL190
FL-410220	220	11.5	12.0	3.2	MSZFL220
FL-410250	250	12.0	12.5	3.4	MSZFL250
FL-410320	320	13.0	13.5	3.8	MSZFL320

STYLE ML MODERATE HEIGHT/LOW PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
ML-410125	125	10.0	9.1	2.9	MSZML125
ML-410170	170	11.0	10.1	3.1	MSZML170
ML-410195	195	11.5	10.6	3.2	MSZML195
ML-410220	220	12.0	11.1	3.4	MSZML220
ML-410285	285	13.0	12.1	3.8	MSZML285

STYLE LL LOW HEIGHT/LOW PROJECTION

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Re-sterilizable Sizer
LL-410135	135	10.5	8.6	3.0	MSZLL135
LL-410180	180	11.5	9.6	3.2	MSZLL180
LL-410210	210	12.0	10.1	3.4	MSZLL210
LL-410240	240	12.5	10.5	3.6	MSZLL240
LL-410300	300	13.5	11.4	4.0	MSZLL300

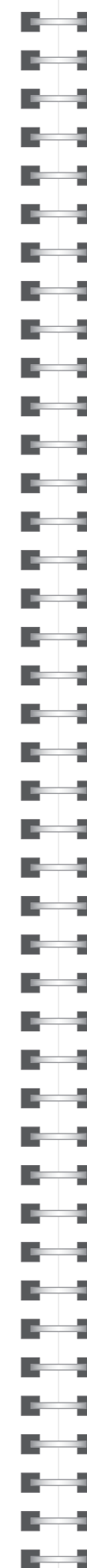

Natrelle[®] **SALINE**



Natrelle[®] **Saline-Filled** Breast Implants— always an option

With *Natrelle*[®] Saline-Filled Breast Implants, you can select an implant that meets your patients' expectations.

***Natrelle*[®] Saline-Filled Breast Implants offer:**

- 98 sizes: 55 smooth and 43 with BIOCELL[®] textured surface³
- Round and shaped options
- 3 projection style options

Natrelle® SALINE-FILLED
SMOOTH ROUND IMPLANTS

STYLE 68HP HIGH PROFILE/ANTERIOR DIAPHRAGM VALVE



Catalog Number	Fill Volume (cc)	Dimensions at Minimum Volume		Dimensions at Maximum Volume		Single-Use Sizer
		Diameter (cm)	Projection (cm)	Diameter (cm)	Projection (cm)	
68HP-160	160-175	9.0	3.9	8.8	4.3	---
68HP-200	200-220	9.6	4.2	9.5	4.6	SZHP68200
68HP-240	240-260	10.2	4.4	10.0	5.0	---
68HP-280	280-300	10.6	4.7	10.5	5.2	SZHP68280
68HP-320	320-340	11.1	4.9	11.0	5.3	SZHP68320
68HP-350	350-380	11.6	4.9	11.4	5.5	SZHP68350
68HP-400	400-430	11.9	5.0	11.7	5.9	SZHP68400
68HP-425	425-455	12.3	5.3	12.1	5.9	SZHP68425
68HP-465	465-505	12.6	5.6	12.5	6.1	SZHP68465
68HP-500	500-540	13.0	5.7	12.8	6.3	SZHP68500
68HP-550	550-590	13.3	5.8	13.1	6.4	SZHP68550
68HP-600	600-640	13.7	6.0	13.5	6.6	---
68HP-650	650-700	14.0	6.1	13.8	6.9	SZHP68650
68HP-700	700-750	14.4	6.2	14.1	7.1	---
68HP-750	750-800	14.6	6.5	14.4	7.2	SZHP68750
68HP-800	800-850	15.0	6.7	14.7	7.2	---

STYLE 68MP MODERATE PROFILE/ANTERIOR DIAPHRAGM VALVE



Catalog Number	Fill Volume (cc)	Dimensions at Minimum Volume		Dimensions at Maximum Volume		Single-Use Sizer
		Diameter (cm)	Projection (cm)	Diameter (cm)	Projection (cm)	
68-120	120-150	9.0	3.0	8.9	3.5	---
68-150	150-180	9.6	3.3	9.5	3.8	---
68-180	180-210	10.2	3.4	10.1	3.9	SZ68180
68-210	210-240	10.6	3.7	10.5	4.1	SZ68210
68-240	240-270	11.1	3.8	11.0	4.2	SZ68240
68-270	270-300	11.6	3.9	11.5	4.3	SZ68270
68-300	300-330	11.9	4.1	11.8	4.4	SZ68300
68-330	330-360	12.3	4.2	12.2	4.5	SZ68330
68-360	360-390	12.7	4.2	12.6	4.5	SZ68360
68-390	390-420	13.0	4.5	12.9	4.6	SZ68390
68-420	420-450	13.4	4.5	13.3	4.6	SZ68420
68-450	450-480	13.7	4.6	13.6	4.7	SZ68450
68-480	480-510	14.1	4.6	14.0	4.7	SZ68480
68-510	510-540	14.4	4.6	14.3	4.7	SZ68510
68-550	550-600	14.6	4.9	14.5	5.1	---
68-600	600-650	15.0	5.0	14.9	5.2	---
68-650	650-700	15.2	5.3	15.1	5.5	SZ68650
68-700	700-750	15.6	5.4	15.5	5.6	---
68-750	750-800	15.9	5.6	15.8	5.8	---
68-800	800-850	16.4	5.6	16.3	5.8	---

Natrelle® SALINE-FILLED
SMOOTH ROUND IMPLANTS


STYLE 68LP LOW PROFILE/ANTERIOR DIAPHRAGM VALVE



Catalog Number	Fill Volume (cc)	Dimensions at Minimum Volume		Dimensions at Maximum Volume		Single-Use Sizer
		Diameter (cm)	Projection (cm)	Diameter (cm)	Projection (cm)	
68LP-125	125-145	9.5	2.6	9.3	3.2	---
68LP-150	150-170	10.0	2.7	9.9	3.2	SZLP68150
68LP-175	175-195	10.6	2.8	10.4	3.3	SZLP68175
68LP-200	200-220	11.0	3.0	10.9	3.4	SZLP68200
68LP-225	225-245	11.4	3.1	11.3	3.4	SZLP68225
68LP-250	250-270	11.9	3.2	11.7	3.5	---
68LP-275	275-295	12.2	3.3	12.1	3.6	SZLP68275
68LP-300	300-320	12.5	3.4	12.4	3.7	SZLP68300
68LP-325	325-345	12.9	3.5	12.7	3.8	SZLP68325
68LP-350	350-370	13.3	3.6	13.2	3.8	SZLP68350
68LP-400	400-420	13.5	3.8	13.4	4.2	SZLP68400
68LP-420	420-440	13.9	3.9	13.8	4.1	SZLP68420
68LP-440	440-460	14.2	3.9	14.1	4.1	---
68LP-480	480-500	14.5	4.0	14.4	4.3	SZLP68480
68LP-525	525-545	14.8	4.2	14.7	4.5	SZLP68525
68LP-550	550-570	15.1	4.3	15.0	4.5	SZLP68550
68LP-600	600-620	15.4	4.4	15.3	4.6	SZLP68600
68LP-640	640-660	15.8	4.5	15.7	4.7	SZLP68640
68LP-680	680-700	16.2	4.5	16.1	4.7	---

STYLE 468 SHAPED

FULL HEIGHT
ANTERIOR DIAPHRAGM VALVE

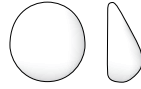

Catalog Number	Fill Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Single-Use Sizer
468-195	195-205	10.0	10.5	4.0	SZ468195
468-230	230-240	10.5	11.0	4.2	SZ468230
468-270	270-285	11.0	11.5	4.3	SZ468270
468-300	300-315	11.5	12.0	4.6	SZ468300
468-350	350-370	12.0	12.5	4.8	SZ468350
468-380	380-400	12.5	13.0	4.9	SZ468380
468-450	450-475	13.0	13.5	5.3	SZ468450
468-495	495-520	13.5	14.0	5.5	SZ468495
468-560	560-590	14.0	14.5	5.7	SZ468560
468-620	620-650	14.5	15.0	5.9	SZ468620

STYLE 363LF SHAPED


LOW HEIGHT/FULL PROJECTION
ANTERIOR DIAPHRAGM VALVE



Catalog Number	Fill Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
363LF-170	170-180	10.4	8.6	4.2
363LF-200	200-210	10.8	9.0	4.5
363LF-230	230-240	11.2	9.4	4.7
363LF-260	260-275	11.8	9.8	4.9
363LF-300	300-315	12.2	10.2	5.1
363LF-330	330-350	12.6	10.6	5.3
363LF-370	370-390	13.0	11.0	5.5
363LF-410	410-430	13.6	11.4	5.6
363LF-450	450-475	14.0	11.8	5.8
363LF-510	510-535	14.4	12.2	6.1
363LF-560	560-585	15.0	12.6	6.2
363LF-620	620-645	15.4	13.0	6.5
363LF-690	690-720	15.8	13.4	6.8

STYLE 168MP ROUND

MODERATE PROFILE
ANTERIOR DIAPHRAGM VALVE

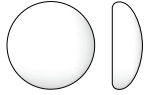

Catalog Number	Fill Volume (cc)	Diameter (cm)	Projection (cm)	Single-Use Sizer
168-120	120-150	9.0	3.0	---
168-150	150-180	9.6	3.3	---
168-180	180-210	10.2	3.4	SZ68180
168-210	210-240	10.6	3.7	SZ68210
168-240	240-270	11.1	3.8	SZ68240
168-270	270-300	11.6	3.9	SZ68270
168-300	300-330	11.9	4.1	SZ68300
168-330	330-360	12.3	4.2	SZ68330
168-360	360-390	12.7	4.2	SZ68360
168-390	390-420	13.0	4.5	SZ68390
168-420	420-450	13.4	4.5	SZ68420
168-450	450-480	13.7	4.6	SZ68450
168-480	480-510	14.1	4.6	SZ68480
168-510	510-540	14.4	4.6	SZ68510
168-550	550-600	14.6	4.9	---
168-600	600-650	15.0	5.0	---
168-650	650-700	15.2	5.3	SZ68650
168-700	700-750	15.6	5.4	---
168-750	750-800	15.9	5.6	---
168-800	800-850	16.4	5.6	---


Natrelle® 133 PLUS TISSUE EXPANDER

Building on the success of the original *Natrelle*® 133— the #1 surgeon-selected tissue expander.^{2,*}

*Based on US market share data through November 2016 (N = 67).

Blue tabs designed for suture accuracy.

Blue orientation line assists in precise placement within the breast pocket.

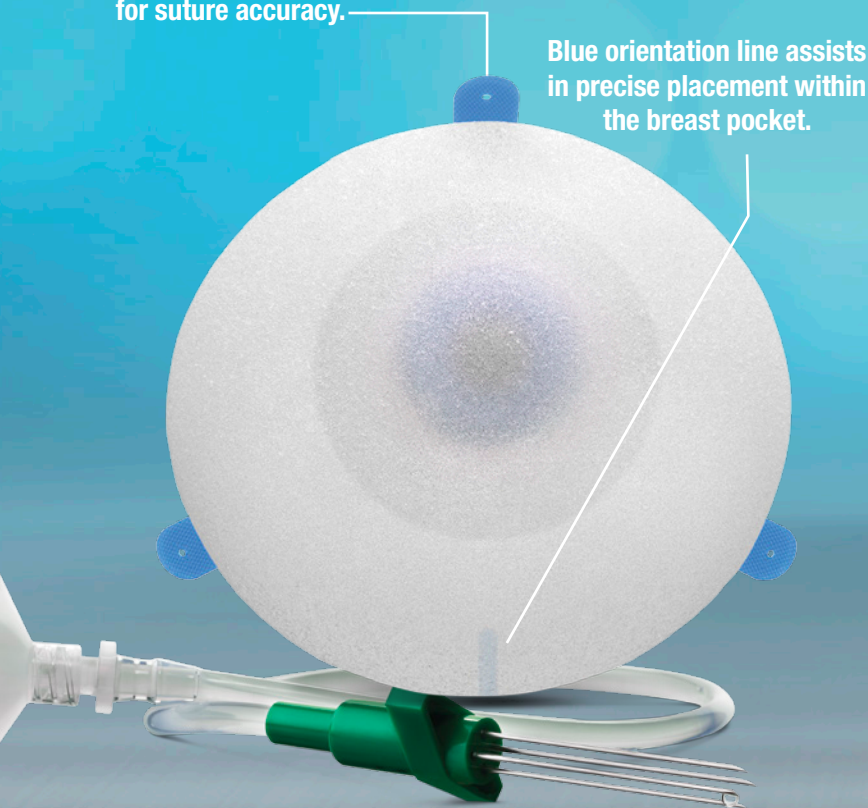

FOURTE™ Expander Fill System

Paired for precision

Natrelle® 133Plus delivers a seamless match between tissue expander and round breast implant.

- 100% match (in base width *and* projection) to *Natrelle* INSPIRA®^{12,†}
- 83% match (in base width *and* projection) to *Natrelle*® 410^{12,†}

Only *Natrelle*® 133Plus Tissue Expanders include the **FOURTE™** Expander Fill System

- Patented 4-needle design fills tissue expanders 4X faster than a standard 21-gauge needle^{9,‡,§}
- May save up to 9 minutes of expander fill time^{9,‡,§}

Natrelle® 133Plus
TISSUE EXPANDERS

STYLE 133P-FX

FULL HEIGHT/
EXTRA PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
133P-FX-11-T	350	11.0	11.5	5.9
133P-FX-12-T	450	12.0	12.5	6.3
133P-FX-13-T	550	13.0	13.5	6.7
133P-FX-14-T	650	14.0	14.5	7.1
133P-FX-15-T	800	15.0	15.5	7.6
133P-FX-16-T	950	15.5	16.0	8.1

STYLE 133P-FV

FULL HEIGHT/
VARIABLE PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
133P-FV-11-T	300	11.0	11.5	5.0
133P-FV-12-T	400	12.0	12.5	5.3
133P-FV-13-T	500	13.0	13.5	5.7
133P-FV-14-T	600	14.0	14.5	6.2
133P-FV-15-T	750	15.0	15.5	6.7
133P-FV-16-T	850	16.0	16.5	6.8

STYLE 133P-MX

MODERATE HEIGHT/
EXTRA PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
133P-MX-11-T	300	11.0	10.0	5.9
133P-MX-12-T	400	12.0	11.0	6.3
133P-MX-13-T	500	13.0	12.0	6.7
133P-MX-14-T	600	14.0	13.0	7.1
133P-MX-15-T	700	15.0	14.0	7.6
133P-MX-16-T	850	15.5	14.5	8.1

STYLE 133P-MV

MODERATE HEIGHT/
VARIABLE PROJECTION




Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
133P-MV-11-T	250	11.0	10.0	4.9
133P-MV-12-T	300	12.0	11.0	5.2
133P-MV-13-T	400	13.0	12.0	5.6
133P-MV-14-T	500	14.0	13.0	6.0
133P-MV-15-T	600	15.0	14.0	6.3
133P-MV-16-T	700	16.0	15.0	6.6

Natrelle® 133Plus
TISSUE EXPANDERS

STYLE 133P-SX

SHORT HEIGHT/
EXTRA PROJECTION

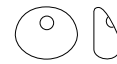

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
133P-SX-11-T	250	11.0	9.0	5.9
133P-SX-12-T	350	12.0	10.0	6.3
133P-SX-13-T	400	13.0	11.0	6.7
133P-SX-14-T	500	14.0	12.0	7.1
133P-SX-15-T	650	15.0	13.0	7.6
133P-SX-16-T	800	15.5	13.5	8.1

STYLE 133P-SV

SHORT HEIGHT/
VARIABLE PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)
133P-SV-11-T	200	11.0	7.6	4.7
133P-SV-12-T	250	12.0	8.4	5.3
133P-SV-13-T	300	13.0	9.1	5.7
133P-SV-14-T	375	14.0	9.7	6.0
133P-SV-15-T	450	15.0	10.5	6.3
133P-SV-16-T	550	16.0	11.3	6.3

STYLE 133P-LV


LOW HEIGHT/
VARIABLE PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm) Tip to Base	Height (cm) Intra-Arc	Projection (cm)
133P-LV-11-T	150	11.0	7.6	7.0	4.7
133P-LV-12-T	200	12.0	8.4	7.5	5.3
133P-LV-13-T	300	13.0	9.1	8.0	5.7
133P-LV-14-T	350	14.0	9.7	8.5	6.0
133P-LV-15-T	400	15.0	10.5	9.0	6.3
133P-LV-16-T	500	16.0	11.3	9.5	6.3

WARNING: The strong rare-earth, permanent magnet contained in the Natrelle® Style 133Plus Series Tissue Expanders is contraindicated where the magnetic field may affect other polarized devices (eg, pacemakers, drug infusion devices, artificial sensing devices, similar-type products, and MRI procedures).

NOTE: All Natrelle® 133Plus Tissue Expanders include the FOURTÉ™ Expander Fill System.

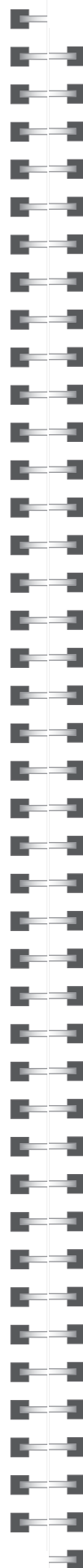

Natrelle® 133
TISSUE EXPANDER

Jessica
Natrelle® 133 Style MV, 400 cc
Natrelle® 410 Style FF, 425 cc
Individual results may vary.


Natrelle® 133
Tissue Expanders

Natrelle® 133

TISSUE EXPANDER

**#1 surgeon-selected tissue expander
with 30 years of clinical experience.**^{2,10,*}

*Based on US market share data through November 2016 (N = 67).

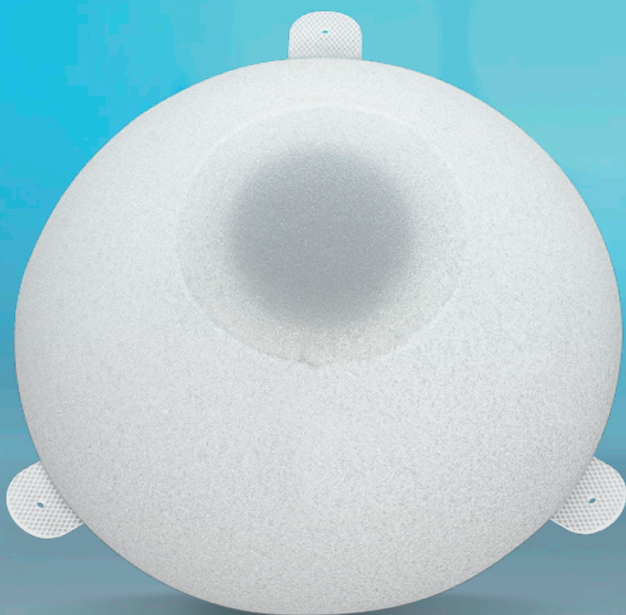

Natrelle® 133 Tissue Expanders

Natrelle® 133 Tissue Expanders offer a precise match with Natrelle® breast implants.


All Natrelle® 133 Tissue Expanders offer:

- 84 options with or without suture tabs³
- BIOCELL® textured surface
- 7 style options

Natrelle® 133
TISSUE EXPANDERS

STYLES 133FX & 133FX-T


FULL HEIGHT/
EXTRA PROJECTION

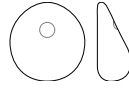

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Catalog Number With Tabs
133FX-11	350	11.0	11.5	5.9	133FX-11-T
133FX-12	450	12.0	12.5	6.3	133FX-12-T
133FX-13	550	13.0	13.5	6.7	133FX-13-T
133FX-14	650	14.0	14.5	7.1	133FX-14-T
133FX-15	800	15.0	15.5	7.6	133FX-15-T
133FX-16	950	15.5	16.0	8.1	133FX-16-T

STYLES 133FV & 133FV-T


FULL HEIGHT/
VARIABLE PROJECTION

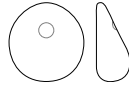

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Catalog Number With Tabs
133FV-11	300	11.0	11.5	5.0	133FV-11-T
133FV-12	400	12.0	12.5	5.3	133FV-12-T
133FV-13	500	13.0	13.5	5.7	133FV-13-T
133FV-14	600	14.0	14.5	6.2	133FV-14-T
133FV-15	750	15.0	15.5	6.7	133FV-15-T
133FV-16	850	16.0	16.5	6.8	133FV-16-T

STYLES 133MX & 133MX-T

MODERATE HEIGHT/
EXTRA PROJECTION

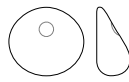

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Catalog Number With Tabs
133MX-11	300	11.0	10.0	5.9	133MX-11-T
133MX-12	400	12.0	11.0	6.3	133MX-12-T
133MX-13	500	13.0	12.0	6.7	133MX-13-T
133MX-14	600	14.0	13.0	7.1	133MX-14-T
133MX-15	700	15.0	14.0	7.6	133MX-15-T
133MX-16	850	15.5	14.5	8.1	133MX-16-T

STYLES 133MV & 133MV-T

MODERATE HEIGHT/
VARIABLE PROJECTION




Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Catalog Number With Tabs
133MV-11	250	11.0	10.0	4.9	133MV-11-T
133MV-12	300	12.0	11.0	5.2	133MV-12-T
133MV-13	400	13.0	12.0	5.6	133MV-13-T
133MV-14	500	14.0	13.0	6.0	133MV-14-T
133MV-15	600	15.0	14.0	6.3	133MV-15-T
133MV-16	700	16.0	15.0	6.6	133MV-16-T

Natrelle® 133
TISSUE EXPANDERS

STYLES 133SX & 133SX-T

SHORT HEIGHT/
EXTRA PROJECTION

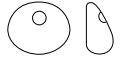

Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Catalog Number With Tabs
133SX-11	250	11.0	9.0	5.9	133SX-11-T
133SX-12	350	12.0	10.0	6.3	133SX-12-T
133SX-13	400	13.0	11.0	6.7	133SX-13-T
133SX-14	500	14.0	12.0	7.1	133SX-14-T
133SX-15	650	15.0	13.0	7.6	133SX-15-T
133SX-16	800	15.5	13.5	8.1	133SX-16-T

STYLES 133SV & 133SV-T

SHORT HEIGHT/
VARIABLE PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm)	Projection (cm)	Catalog Number With Tabs
133SV-11	200	11.0	7.6	4.7	133SV-11-T
133SV-12	250	12.0	8.4	5.3	133SV-12-T
133SV-13	300	13.0	9.1	5.7	133SV-13-T
133SV-14	375	14.0	9.7	6.0	133SV-14-T
133SV-15	450	15.0	10.5	6.3	133SV-15-T
133SV-16	550	16.0	11.3	6.3	133SV-16-T


STYLES 133LV & 133LV-T

LOW HEIGHT/
VARIABLE PROJECTION



Catalog Number	Volume (cc)	Width (cm)	Height (cm) Tip to Base	Height (cm) Intra-Arc	Projection (cm)	Catalog Number With Tabs
133LV-11	150	11.0	7.6	7.0	4.7	133LV-11-T
133LV-12	200	12.0	8.4	7.5	5.3	133LV-12-T
133LV-13	300	13.0	9.1	8.0	5.7	133LV-13-T
133LV-14	350	14.0	9.7	8.5	6.0	133LV-14-T
133LV-15	400	15.0	10.5	9.0	6.3	133LV-15-T
133LV-16	500	16.0	11.3	9.5	6.3	133LV-16-T

WARNING: The strong rare-earth, permanent magnet contained in the Natrelle® Style 133 Series Tissue Expanders is contraindicated where the magnetic field may affect other polarized devices (eg, pacemakers, drug infusion devices, artificial sensing devices, similar-type products, and MRI procedures).

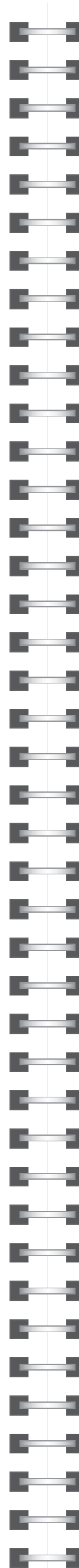

Natrelle[®] tissue expander accessories

Sizer Templates

Description	Catalog Number
Plastic measuring templates for Styles 133Plus and 133 Tissue Expanders, nonsterile:	
Styles 133P-FX and 133FX Sizer Template Set	30-00038
Styles 133P-FV and 133FV Sizer Template Set	30-00018
Styles 133P-MX and 133MX Sizer Template Set	30-00039
Styles 133P-MV and 133MV Sizer Template Set	30-00017
Styles 133P-SX and 133SX Sizer Template Set	30-00040
Styles 133P-SV and 133SV Sizer Template Set	30-00035
Styles 133P-LV and 133LV Sizer Template Set	30-00016

MAGNA-FINDER[®] Xact Port Finders

Description	Catalog Number
Xact2 External locating device with dermal indicator for use with Style 133Plus Tissue Expanders for locating MAGNA-SITE [®] injection sites (provided nonsterile)	30-00034
Xact External locating device for use with Style 133 Tissue Expanders for locating MAGNA-SITE [®] injection sites (provided nonsterile)	30-00032

FOURTE[™] Expander Fill System

Description	Catalog Number
Four-pronged 21-gauge, 1¾" butterfly needles, tubing, luer adapter	F-4444

21-G Needle Infusion Set

Description	Catalog Number
21-gauge, 1¾" butterfly needle, tubing, luer adapter	30-00012

Universal Fill Kit

Description	Catalog Number
60-cc syringe, 122-cm transfer set with piercing device, 25-cm extension tube, 2-way check valve	30-00033

The above accessories are not eligible for return for credit.

Natrelle[®] breast implant accessories

Diaphragm Valve Fill Tube and Reflux Valve

Description	Catalog Number
Diaphragm valve fill tube and reflux valve	25-00017

BIOCELL[®] Delivery Assistance Sleeve

Description	Catalog Number
Clear, sterile sleeve	27-000001

The above accessories are not eligible for return for credit.

Warranty program brochure



ConfidencePlus[®] warranty program

Get more coverage and confidence with *Natrelle*[®] Gel implants than with any other breast implant warranty—**ONLY** with *Natrelle*[®] *ConfidencePlus*[®].

- The **ONLY** warranty with a 10-year FREE breast implant replacement in the event of capsular contracture (Baker Grade III/IV)
- The **ONLY** warranty to cover primary AND revision surgeries in the event of capsular contracture (Baker Grade III/IV)
- The **ONLY** warranty that allows you to choose ANY style and size of *Natrelle*[®] Gel breast implant for replacement in the event of capsular contracture (Baker Grade III/IV) or rupture
- Lifetime product replacement for all gel implants in the event of rupture
- Financial assistance up to \$3500 for 10 years in the event of implant rupture

For more information about the *ConfidencePlus*[®] warranty program, please visit Natrelle.com/warranty.

Additional information

Committed to advancing the field of breast aesthetics.

We have advanced the medical aesthetics industry through scientific discovery, market innovation, and partnering with the medical communities that we serve.

Building on 50 years of scientific research and development, we are taking our passion for innovation and applying it to the consumer realm. That's why we've created the *Natrelle*® Portfolio—specially designed to give you and your patients more options to help you achieve the desired results.



CAUTION: United States Federal Law restricts these devices to sale by, or on the order of, a licensed physician.

PRODUCT SUPPLIED STERILE: Please see Directions for Use at www.allergan.com/labeling/usa.htm for additional product information including warnings, precautions, adverse reactions, and instructions for use. *Natrelle*® Styles 133Plus and 133 Tissue Expander Sizing Templates and MAGNA-FINDER® Xact2 and Xact are supplied nonsterile.

DIMENSIONS: Dimensions listed are typical and are measured with the implant placed on a flat surface. Not every implant will conform to the dimensions given. Slight variations may occur.

RECOMMENDED FILL VOLUMES: Following recommended fill volumes for saline-filled breast implants can decrease the possibility of shell wrinkling and crease-fold failure.¹³

Shaping the future with you

Natrelle® is the #1 plastic surgery portfolio.^{1,2,*†}

Natrelle INSPIRA® Collection

#1 selected round gel implant collection in the US.^{1,*}

Natrelle® 410

#1 selected anatomically shaped implant in the US.^{2,†}

Natrelle® 133 Series

#1 selected tissue expander in the US.^{2,†}

Natrelle INSPIRA® Cohesive

Most cohesive round gel breast implant in the US.¹¹

*Based on surgeon survey data, December 2016 (N = 100).

†Based on US market share data through November 2016 (N = 67).

To order breast aesthetics products,
please call Customer Care at 1-800-766-0171.

ART. SCIENCE. *Natrelle*®

References:

1. Data on file, Allergan, December 2016; Plastic Surgery Monthly Tracker: Breast Implants & Tissue Expanders.
2. Data on file, Allergan, November 2016; Breast Implant + Tissue Expander Share Tracker.
3. More of Everything: *Natrelle*® Product Catalog. Irvine, CA: Allergan; 2015.
4. Our products: breast implants. Mentor Worldwide LLC website. <http://www.mentorwwllc.com/global-us/Breast.aspx>. Updated April 2016. Accessed February 14, 2017.
5. *Sientra*® HSC & HSC + Breast Implant Product Catalog. Santa Barbara, CA: Sientra, Inc.; 2015. <http://sientra.com/Content/pdfs/MDC-0177%20R2%20HSC%20HSC%20Breast%20Implants%20Product%20Catalog.pdf>. Accessed February 14, 2017.
6. Data on file, Allergan, January 2017; Allergan Plastic Surgery Order Form.
7. *Natrelle*® Silicone-Filled Breast Implants and *Natrelle* INSPIRA® Breast Implants: Smooth & BIOCELL® Texture Directions for Use, 2016.
8. *Natrelle*® 410 Highly Cohesive Anatomically Shaped Silicone-Filled Breast Implants Directions for Use, 2014.
9. Data on file, Allergan, January 4, 2016; Protocol MM-1225-FR.
10. Data on file, Allergan, July 14, 1986; FDA Section 510(k) marketing approval letter.
11. Data on file, Allergan, February 8, 2017; Study Report MD16075-DV2.
12. Data on file, Allergan, January 6, 2017; Study Report MD16076-DV.
13. *Natrelle*® Saline-Filled Breast Implants: Smooth & BIOCELL® Texture Directions for Use, 2014.



© 2017 Allergan. All rights reserved. All trademarks are the property of their respective owners.
Natrelle.com NAT64504_v2 03/17 170324